

# Operating table RiEye Mk2S Accessories



## RiEye Mk2S

Rini's operating tables in the "RiEye"-series are popular for microsurgery in both in stationary use and in "roll-in roll-out" operation. They are easy to handle and have many useful accessories. The stable construction is ergonomically designed to provide a high level of comfort, both for the patient and medical personnel

**RiEye Mk2S comes in three different base models;**



R5 model enables the patient to sit down or stand up from the front of the operating table



R6 is popular for cataract surgery, IVT eye injections, plastic surgery and similar microsurgical procedures



R7 operating table is suited for hospitals and clinics where surgeries under full anaesthesia are performed

### GENERAL

Two column electrical system with adaptive chassis

Weight of base unit	95kg /209lb	Model R5 length	1700-1940mm/67-76in
Lifting capacity patient	300kg/661lb	Model R6 length	1790-2030mm/70-80in
Power system	Battery	Model R7 length	1960-2200mm/77-87in
External charger	110-240V	Footprint wheelbase length x width	800 x 600mm/31-23in
Height	500-900mm/19-35in	Cushion width	580mm/23in
Height with reverse Trendelenburg	550-950mm/21-37in	Cushion width extra wide	680mm/27in

### COLORS



### RiEye Mk2S accessories

There are a variety of accessories that have been developed in close cooperation with hospitals and clinics. These enable optimization of different treatment methods and provide an ergonomic and comfortable working environment. Last but not least, they also contribute to provide maximum safety and comfort for the patient.

# HEADRESTS AND HEADREST ACCESSORIES

## Doublelink Bowl headrest included in the RiEye

Keeps the head comfortably fixed and is especially suitable for eye surgery



**140-00270-00 Headrest X-tray**  
Large tray with various pillows



**140-00270-01 Pillow for X-tray black**  
Soft memory foam, allowing flexible head position



**140-00270-02 Pillow for X-Tray grey**  
Moulded, cupped shaped soft pillow keeps the head comfortably fixed



**140-00252-05 Headrest FaceDown/FaceUp**  
Tray that allows the patient to lie on the stomach or the back




**140-00270-06 Pillow for FaceDown/FaceUp**  
Moulded, soft pillow both for facedown or faceup surgeries



**140-00270-03 Headrest S-tray**  
Small tray with various pillows



**140-00270-04 Pillow for S-tray**  
Soft memory foam, good for when the head is tilted to the side

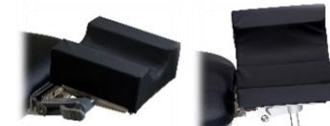

**140-00270-14 Pillow for S-tray, U-shaped**  
Soft memory foam, with firmer sides for a comfortable fixed head position



**140-00270-07 Headrest One grip**  
Tiltable headrest with seamless adjustment in angle, lateral tilt. With a quicklock to change between different headrests



**140-00270-08 Luigi headrest with pillows, for One grip headrest**  
With multiple replaceable cushions in different thickness/firmness



**140-00270-09 S-tray, for onegrip**  
Small tray for various pillows. The pillows are attached with velcro



**140-00229-01 electrical headrest for Doublelink**  
Gives the surgeon the possibility to adjust the headrests angle by the hand- or foot control. For both standard and optional headrests. Can't be retrofitted



**140-00242-01 Wrist support Bowl / S-Tray for Doublelink**

**140-00242-02 Wrist support X-Tray for Doublelink**

Provides relief for the surgeon near the patient's head. Foldable and adjustable in height and length. Can easily be mounted and removed if required



**140-00235-02 Accessory rail for X-tray**

For when accessories will be mounted on the headrest, such as wrist support



**140-03420-00 Bowl headrest (standard) with foundation to make headrest dockable**

Complete headrest foundation with Easy fix headrest (or the ordered headrest), when multiple headrests can be used on the same table

## HAND CONTROL HOLDER



### **140-00238-01 Flexible, bracket included**

Rotatable on flexible bar and is mounted with clamping bracket on the wheel frame



### **140-00238-02 Fixed, bracket included**

Fixed to bar and can be rotated and adjusted vertically



### **140-00238-03 Fixed for IV-pole, bracket included**

Fixed on IV pole and can be rotated and adjusted vertically

## FOOT CONTROL

With four functions: Height, backrest tilt, rotation and leg support tilt or electrical headrest (if electrical headrest option chosen). Allows the surgeon to make fine adjustments during surgery.



### **140-00230-01 Foot control with two joysticks movable tray**

Placed on extendable tray or can be placed on the floor



### **140-00230-02 Foot control with two joysticks fixed tray, can't be combined with electric headrest**

Placed on a fixed tray or on the floor. Tray is higher up than the movable tray to keep access when other pedals coming close



## PERSONAL COMFORT



### **140-00261-01 Auto access to memory location**

Makes the table go automatically to the desired position with one touch of button 1-4, Trend or Level-0.



### **140-00266-01 Electrical brake and direction control**

Brake and direction control is operated by the hand control instead of the brake pedal, increases space next to the chassis. Can't be retrofitted



### **140-00352-00 Reverse Trendelenburg max 12°**

The whole operating table can be tilted forward up to 12 degrees. Seat, leg and back can be individual adjusted while in Reverse Trendelenburg position. Can't be retrofitted

## GENERAL ACCESSORIES



**140-00236-02 Drape support, formable**  
Holds up the surgical drape from the patient.  
Flexible rod that easily bends to the correct shape



**140-00236-03 Drape support, formable with air tube connection**  
Holds up the surgical drape from the patient.  
Flexible rod that easily bends to the correct shape. The nozzle can be rotated to the desired exhaust direction.



**140-00235-00 Attachment rail for accessories on the seat 250mm/10in**  
Convenient if accessories need to be mounted by the seat. The attachment rail can be mounted on the right and / or left side of the seat



**140-00235-01 Attachment rail for accessories on the seat 250mm/10in if Side rail is mounted**  
Convenient if accessories need to be mounted by the seat. The attachment rail can be mounted on the right and / or left side of the seat



**140-00163-00 Anaesthesia armrest, clamp and arm strap included**  
Angle and inclination can be adjusted with a one-hand grip. Can be placed on the left or right side. The arm cushion can be removed for cleaning



**140-00239-01 IV pole, bracket included**  
The IV pole is connected directly to the operating table and can be mounted on or off when needed. Can be placed on the right or left side of the operating table and can be raised and lowered with an ergonomic knob



**140-00237-02 Battery lead (Standard)**

**140-00237-04 Battery Lithium iron.**

Twice the number of operating cycles and half the weight of the standard, and a longer lifespan



**140-00251-00 Leg cushion cover R6 R7**

**140-00251-01 Leg cushion cover R5**

Protects against dirt from shoes. Extends the lifespan of the leg cushion



**140-00351-00 Bracket for oxygen tank**

Provides availability of oxygen when the operating table is mobile

Length of oxygen tank up to 800 mm/31in. Diameter of oxygen tank up to 120mm/5in. Weight of oxygen tank max 10 kg/22lb

## PATIENT COMFORT



**140-00350-01 Width extension back rest left**

**140-00350-00 Width extension back rest right**

Improves patient comfort and supports their arms to be on the operating table and gives a more secure feeling



**140-00232-00 Length extension back rest**

**120mm/5in**

Length extension increases comfort for taller patients



**140-00265-00 Leg support circular**

**Diameter 150mm/6in**

Increased comfort for patients if their knees must come up higher



**140-00265-01 Leg support half circular**

**200mm/8in high**

## TRANSPORT OF THE TABLE



### **140-00259-00 Driving handles backrest foldable R/L**

The driving handles at the backrest fits all models of the operating table and are foldable



### **140-00256-00 Driving handles leg foldable R/L R6 and R7**

Facilitates patient transport. The driving handles are foldable



### **140-00256-01 Driving handles leg R/L R5**

For when the operating table is used as a mobile unit, or it generally needs to be moved

## PATIENT SAFETY



### **140-00240-00 Side rail collapsible right side**

### **140-00240-01 Side rail collapsible left side**

Increased patient safety, especially when moving the operating table. Foldability enables the patient to sit down or stand up from the side of the operating table



### **140-00320-00 Armrest for Side rail collapsible**

Improves comfort for the patient, mounted on the collapsible side rail



### **140-00240-02 Side rail right or left side 600mm/24in, 140-00235-00 included**

Increased patient safety, especially when moving the operating table. Can be mounted on the right and / or left side. Large fixed side rail that is adjustable in height



### **140-00263-00 Belt waist**

Increases patient safety. Easy to install and remove. Wide belts with secure Velcro fastening




### **140-00263-01 Belt arm**

700mm/28in



### **140-00263-02 Belt leg**

2100mm/83in

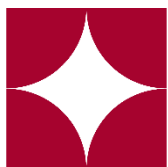

rini

Rini Ergoteknik AB  
Fågelsångsvägen 3B  
SE-186 42 Vallentuna  
Sweden  
Phone +46 8 594 77170  
Email [info@rini.se](mailto:info@rini.se)  
[www.rinicompany.com](http://www.rinicompany.com)



MEDICONTUR  
— AVI —

Medicontur AVI

España y Portugal

Tel.: +34 91 676 47 97  
Email: [info@medicontur.es](mailto:info@medicontur.es)  
[www.medicontur.es](http://www.medicontur.es)